



TheraLux 10K Lamp



Keep Health Always Around You

Thank You.

Thank you for purchasing the TheraLux 10K Lamp
We are dedicated to providing our customers with quality products for building a health living. Should you have any questions or concerns about using your new product, feel free to reach out to our helpful customer support team We hope you enjoy your new product!

IMPORTANT SAFETY INSTRUCTIONS

Read and follow all instructions in this user manual and heed all warnings

WARNING

To reduce the risk of burns,fire,electric shock,or injury to persons:

- This appliance is for indoor use only.
- Water and electricity are a dangerous combination. Do not use this appliance in wet surroundings(e.g. in the bathroom or near a shower or swimming pool).
- Make sure the adapter does not get wet.
- Do not let water run into the appliance or spill water onto the appliance.
- Do not immerse the appliance in water or rinse under the tap.
- The adapter contains a transformer, which converts an unsafe 100-240V AC to a safe 24V DC low voltage. Do not cut off the adapter to replace it with another plug, as this causes a hazardous situation.
- Do not use the adapter in or near a power outlet that contains an electric air freshener to prevent damage to the adapter.
- Protect the power cord from being walked on or pinched, particularly at plugs, outlet and point where they exit from lamp.
- Do not impair the cooling of the appliance by covering the appliance with items such as blankets, curtains, clothes, papers, etc.

CAUTION


- Do not drop the appliance on the floor, hit it hard or expose it to other heavy shocks.
- Place the appliance on a stable, level and non-slippery surface.
- This appliance is only intended for household use, including similar use in hotels.
- Check if the voltage indicated on the appliance corresponds to the local voltage before you connect the appliance.
- Unplug this appliance during Lightning storms or when unused for long periods of time.

GENERAL

- Consult your doctor before you start using the appliance if you have suffered from or are suffering from depression.


-1-

CONTENTS

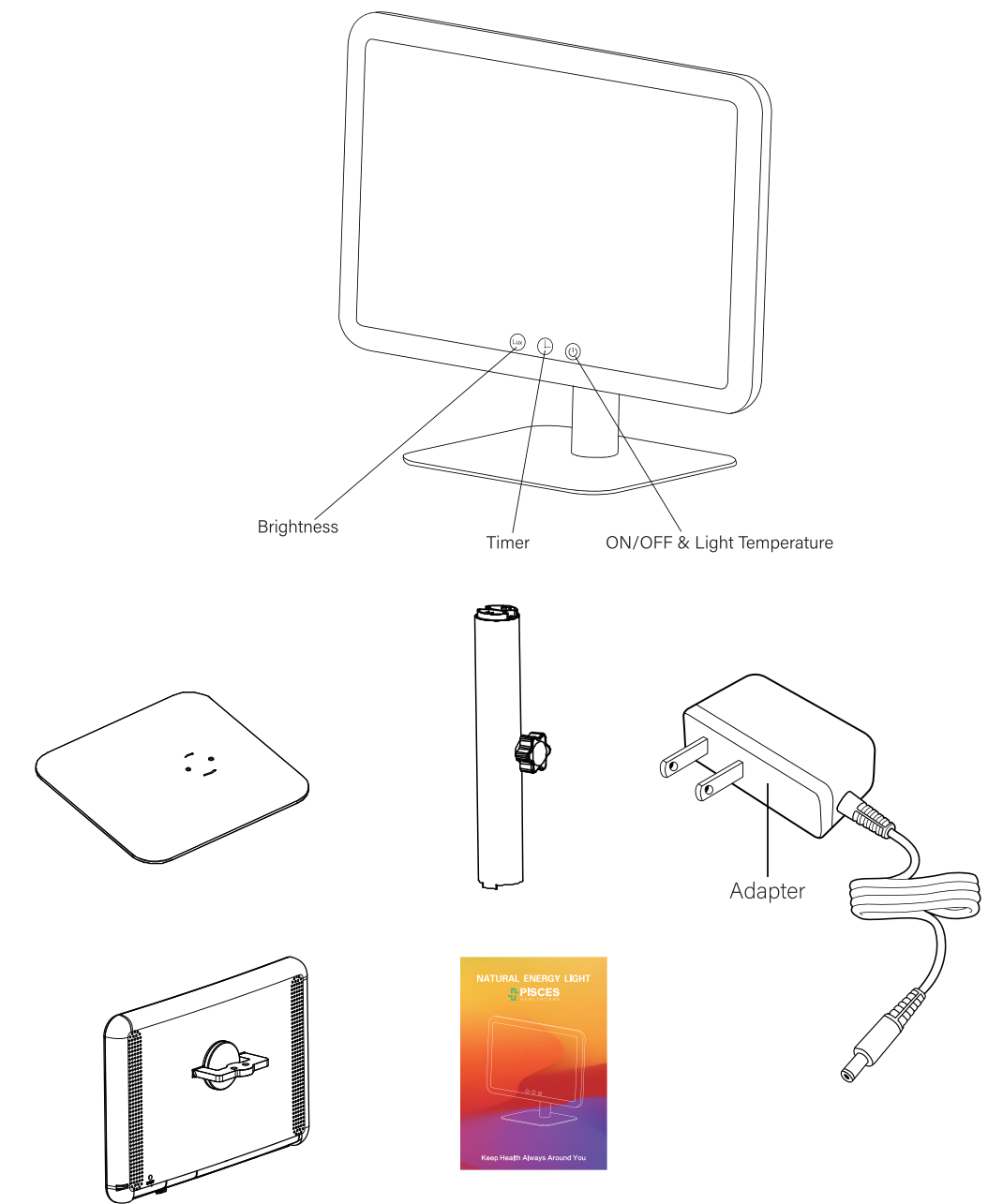

-2-

ASSEMBLE

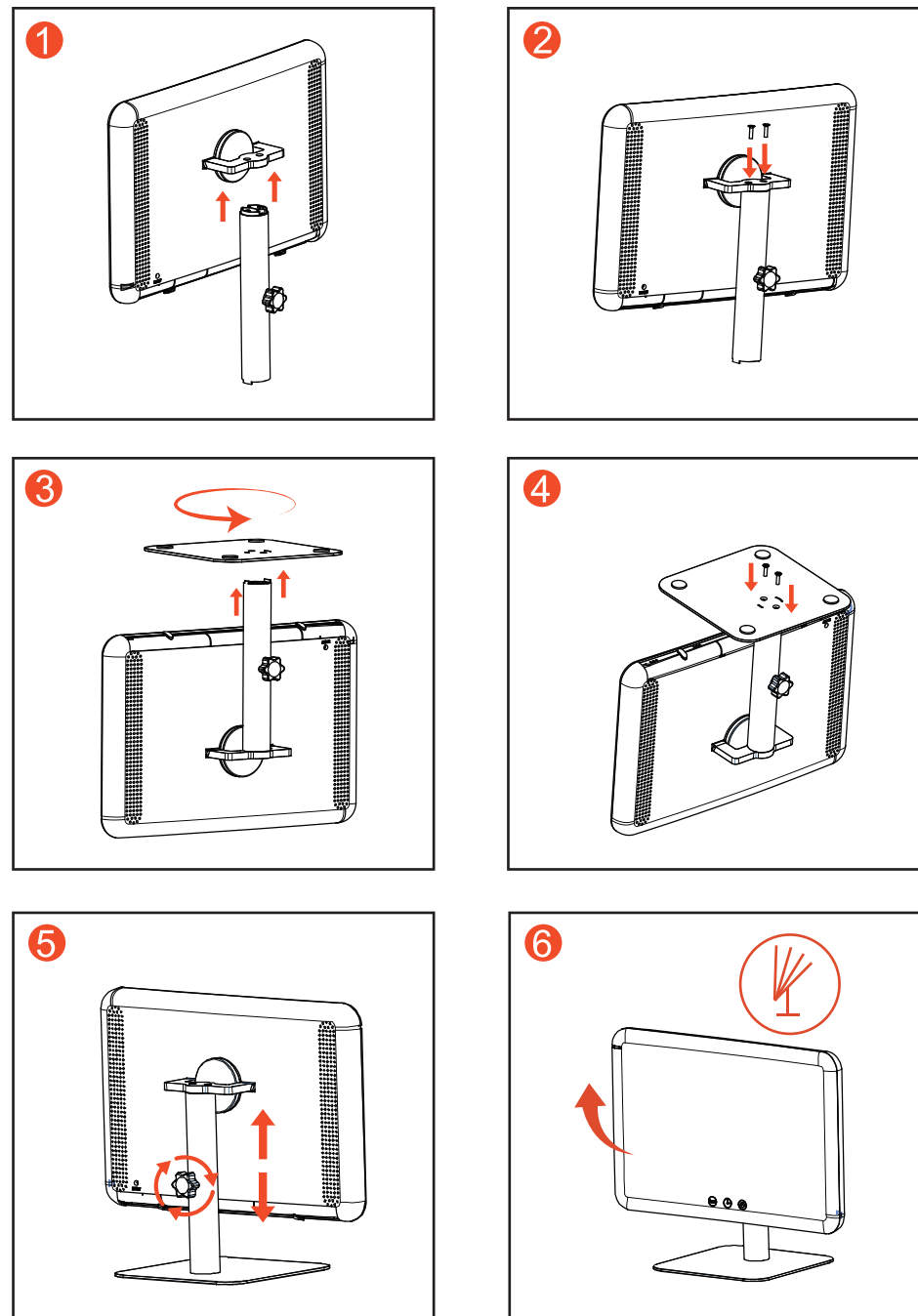

-3-

FUNCTIONS

- Short press to turn on the light
- Long press the color temperature key to adjust different color temperature modes
- Press the timer button when the light on, 15min-30min-60min-90min
When the light on it means timer set ok
And press the button again to cancel the timer
- Press the brightness button to adjust the brightness of the light, the ratio is 20%-40%-60%-80%-100%

-4-

SPECIFICATIONS

Model:	TheraLux 10K Lamp
Light color:	3000-6500K
Timer:	15-30-60-90min
Rated output power adapter:	24W
Light output:	>10000Lux
Rated input voltage adapter:	100-240VAC



-5-